Automated Time-Lapse Microscopy — Inženernoe reshenie + budget

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Компонент	Модель / источник	Цена
Motorized XY stage adapter	Arduino CNC-shield + stepper motors + linear rails (custom) or OpenFlexure-style DIY kit	\$500
Motorized focus drive	NEMA-17 stepper + belt drive on Zeiss coarse knob (3D-printed mount)	\$200
Digital camera	FLIR Blackfly S USB3 BFS-U3-63S4M-C (mono, 2448×2048, scientific)	\$1,200
Live cell chamber	OkoLab bold-line стандарт (\$1,500) OR DIY acrylic chamber + Peltier heater + CO ₂ solenoid (\$400)	\$400 (DIY) or \$1,500 (pro)
CO₂ + temperature controller	OkoLab H301 or Tokai Hit INU compact (\$800) OR DIY with CO ₂ tank + Arduino + DHT22 sensor + solenoid (\$250)	\$250 (DIY)
LED fluorescence illuminator	SOLA SE II Lumencor (\$1,500 used) or ThorLabs M470L4 + M565L3 + filter cube (\$500)	\$500 (LED stand-alone)
Computer	Refurbished Dell OptiPlex i5 16GB RAM 512GB SSD (\$300)	\$300
Acquisition software	Micro-Manager 2.0 (free, open source) + PyMMCore-Plus Python bindings	\$0
Data processing	Fiji (free) + custom Python pipeline	\$0
UPS	APC Back-UPS Pro 1500VA (\$200)	\$200
Remote monitoring	VPN (WireGuard free) +	\$0

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Компонент	Модель	Цена
Prior 3rd-party motorized XY stage	Prior ProScan III H117 с controller	\$3,500 (used eBay)
Prior motorized focus (Z)	Prior Z-Drive motor + adapter	\$1,500 (used)
Digital camera	Hamamatsu ORCA-Fusion BT (\$6k new, \$2k used)	\$2,000 (used)
Live cell chamber + CO₂	OkoLab H301-K-Frame bold-line complete	\$2,500 (used)
LED fluorescence	Lumencor SOLA SE II	\$1,500 (used)
Computer	Dell Precision 3650 Xeon, 32GB RAM, 1TB NVMe	\$800
Software	Micro-Manager (free) or µManager Studio	\$0
UPS + NAS + miscellaneous		\$900
Подитог Вариант В		\$12,700

Плюсы: Industrial-grade reliability. 2-3 days setup. No calibration issues.

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Philosophy: Buy ready-to-use automated microscopy system, ignore Zeiss IM 35.

- **Nikon Eclipse Ti2-E** motorized inverted + Perfect Focus System (PFS) + DS-Qi2 camera + NIS-Elements HCA-module + OkoLab environmental
- Used market price: \$25,000–50,000 depending on age and configuration

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Новая budget line в LOI Phase A

Раздел “Microscope Upgrade & Automation”: \$4,500 (было \$2,000)

Sub-item	Cost
Motorized XY stage (DIY kit / Arduino CNC)	\$500
Motorized Z focus drive	\$200
FLIR Blackfly S digital camera	\$1,200
LED fluorescence illuminator (ThorLabs M470L4+M565L3)	\$500
CO ₂ + environmental DIY controller	\$250
DIY live cell chamber (acrylic + heater)	\$400
Computer (refurbished)	\$300
UPS 1500VA	\$200
NAS backup system	\$500
3D-printed mounts + cables + misc	\$150
Setup contingency	\$300
Total	\$4,500

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Timeline for setup

Week	Task
1	Ordering components (Arduino, FLIR, chamber, LED, NAS)
2	Delivery + initial assembly
3	Motorized stage integration + Z-focus motor mounting
4	Camera + Micro-Manager installation + calibration
5	Environmental chamber + CO ₂ flow testing
6	LED fluorescence setup + filter alignment
7	Full system integration test (dry run without cells)
8	First cell experiment (BJ-hTERT validation)
9-24	Phase A data collection (experiments Aim A.1, A.2, A.3)

Setup occupies Months 1-2 of Phase A budget. Experiments begin Month 3. Data analysis through Month 6.

Technical specifications

Imaging parameters

Parameter	Value
Resolution	2448×2048 px (FLIR BFS-U3-63S4M)
Bit depth	12-bit → 16-bit stored
Frame rate	74 FPS max (not limiting for time-lapse)
Exposure	50-500 ms (brightfield), 500-3000 ms (fluorescence)
Z-stack	5-11 slices × 2 μm spacing (optional for 3D imaging)
Channels	Brightfield + FITC (GT335/polyE) + TRITC (Ninein)
Time-lapse interval	30 min (day/night)
Data per session	~5 GB/day (if 10 positions × 3 channels × 48 timepoints)
Total data 6 mo	~900 GB (easily fits NAS 4TB)

Stage automation

- **XY range:** Ensure coverage of 35mm dish ($\sim 40\text{mm} \times 40\text{mm}$) or 6-well plate
- **XY accuracy:** $\pm 5 \text{ }\mu\text{m}$ positioning (achievable with Arduino CNC-style steppers + linear rails; Prior stages are $\pm 0.1 \text{ }\mu\text{m}$ but more expensive)
- **Z range:** 10mm travel, $\pm 0.5 \text{ }\mu\text{m}$ accuracy (adequate for autofocus)
- **Autofocus:** Software-based, contrast detection on DAPI/Hoechst channel OR hardware laser-based (more expensive, Перспектива в upgrade)

Environmental control

- **Temperature:** $37^\circ\text{C} \pm 0.3^\circ\text{C}$ (Peltier + PID controller)
 - **CO₂:** $5\% \pm 0.2\%$ (mass flow controller + CO₂ tank; one tank lasts ~ 2 weeks)
 - **Humidity:** 95% (passive evaporation from humidifier insert)
-

Remote monitoring setup

1. **VPN** (WireGuard) into lab network from Jaba's home
 2. **Micro-Manager live display** via VNC/RustDesk
 3. **Camera preview stream** via MJPEG over HTTP (for quick peek)
 4. **SMS alerts** via Twilio API ($\$0.01/\text{msg}$, $\sim \$1/\text{month}$) when:
 - o Temperature deviates $> 0.5^\circ\text{C}$
 - o CO₂ deviates $> 0.5\%$
 - o Camera disconnects
 - o UPS engaged (power failure)
 5. **Daily email summary** — automated Python script sends overview image + environmental log at 8am
-

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Camera + optics (\$1,700 total)

- FLIR Blackfly S USB3 BFS-U3-63S4M-C monochrome scientific camera (\$1,200, new)
- C-mount to Zeiss IM 35 adapter (\$80)
- ThorLabs M470L4 blue LED (470 nm, for FITC excitation) (\$180)
- ThorLabs M565L3 green LED (565 nm, for TRITC excitation) (\$180)
- ThorLabs LEDD1B LED driver (\$60)

Live cell environmental (\$1,050 total)

- Acrylic chamber (custom-cut 4mm thick, ~200mm × 200mm × 60mm internal, \$120)
- Glass cover slip for bottom (\$20)
- CO₂ tank 9kg (refillable, \$60 Tbilisi medical supply)

- CO₂ regulator + flow meter (\$80)
- Peltier element TEC1-12706 + heatsink + fan (\$25)
- PID controller (Inkbird ITC-100 or equivalent) (\$35)
- Humidifier insert (plastic tray with foam) (\$10)
- OKolab compact CO₂ incubator mini chamber (used market, optional upgrade, \$700)

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- Уменьшить "General Consumables" с \$12,000 на \$9,500 (-\$2,500)
- Увеличить "Microscope Upgrade" с \$2,000 на \$4,500 (+\$2,500)
- Переименовать line в "**Microscope Automation & Upgrade**"
- Добавить footnote: "Includes motorized XY + Z stage for round-the-clock imaging, FLIR digital camera, LED fluorescence illumination, environmental chamber (37°C + 5% CO₂), NAS backup, and UPS. Full specification in supplementary materials."

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⚠ Key risks & mitigations

1. **DIY stage calibration drift** → use linear rails с encoders (add ~\$100); alignment every 2 weeks

2. **CO₂ tank run-out mid-experiment** → install low-pressure sensor + SMS alert; keep 2 tanks rotating
 3. **Power outage > UPS runtime** → UPS gives 2h; manually restart if needed; alert sent immediately
 4. **Contamination** (long-term cultures) → weekly media change, sterile technique, monthly bleach decontamination of chamber
 5. **Camera failure** → budget \$500 contingency for replacement (not in main budget; covered by \$300 setup contingency + technician hours)
 6. **Software crash** → Micro-Manager + watchdog script: if no new frame in 2h, restart acquisition automatically
 7. **Remote access compromised** → WireGuard with ed25519 keys; no public SSH; lab IP whitelist
-


Claude Code /overnight — AI-operated microscopy

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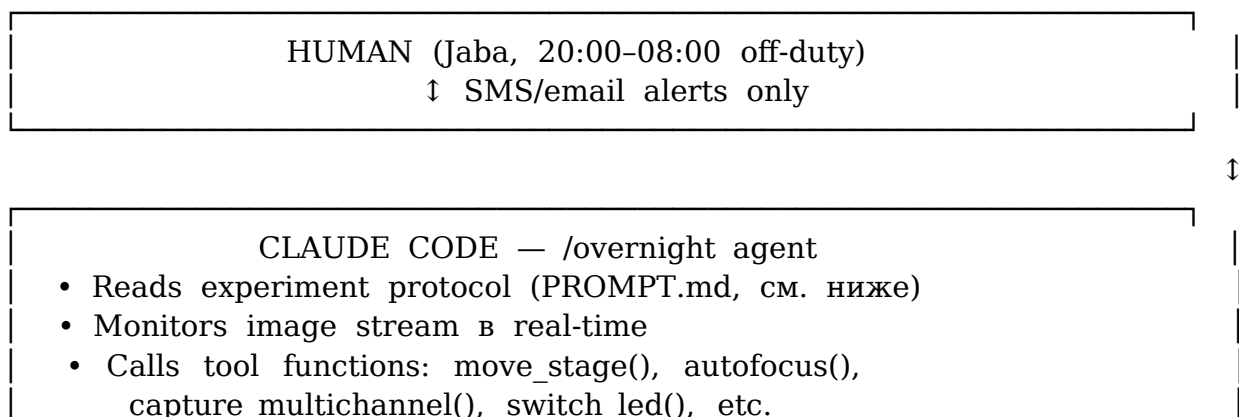- Decision loop: каждые 30 мин inspect last batch, выбрать next action (stay/move/rotate/alarm)
- Journalled decisions с rationale
- On blocker: pauses + SMS to experimenter

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```
move_stage_xy(x_mm, y_mm); move_stage_relative(dx, dy)
autofocus(mode='contrast'|'laplacian')
capture_single(channel, exposure_ms); capture_zstack(ch, z_range, step)
capture_multichannel(channels)
switch_led(channel, intensity_pct)
```

#

image

analysis

```
detect_cells(image); detect_centrioles(image)
measure_polyglu_intensity(cell_roi, channel)
segment_nuclei(image); count_dead_cells(image)
```

#

environmental

```
get_temperature_c(); get_co2_percent(); get_humidity_percent()
set_target_temp(c); open_co2_valve(duration_s)
```

#

alerting

/

journaling

```
send_sms_to_human(message, priority='info'|'warn'|'crit')
log_decision(action, rationale, observed)
save_checkpoint(state); halt_experiment(reason) # requires human unlock
```

Cost

Claude Code subscription: **\$20/мес Pro tier** × 6 мес = **\$120 total** — negligible vs \$4,500 hardware.

PROMPT-based experimentor signaling

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~/Documents/Experiments/CDATA_PhaseA_Aim1/PROMPT.md:

Aim A.1 — α vs β discrimination

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Disentangle division-dependent (α) from time-dependent (β) polyGlu accumulation in BJ-hTERT fibroblasts by comparing proliferating vs contact-inhibited conditions.

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- ****Proliferating wells**** (A1-A6): sub-confluent, passaged every 3 days
- ****Contact-inhibited wells**** (B1-B6): seeded confluent, media changed 2× per week
- ****Time points for imaging****: days 0, 7, 14, 21
- ****Channels****: Brightfield, FITC (GT335 polyGlu), TRITC (Ninein mother centriole)

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- polyGlu intensity at mother centriole (FITC co-localized with TRITC)
- cell density (brightfield segmentation)
- morphology changes (senescence markers: flattened cells, vacuoles)

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- ***Experiment compromised***: если cell death >40% in any well → stop imaging, SMS Jaba, save last data
- ***Environment failure***: если temp <35°C, CO₂ <4% или >6%, humidity <80% — pause, SMS Jaba
- ***Hardware fault***: если camera disconnects, stage jams, UPS engaged — halt, SMS Jaba

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- ****Daily 8 AM email****: summary of last 24 hrs + 4 sample images
- ****Mid-day SMS (14:00)****: short status (OK / warning / critical)
- ****Cycle log (every 30 min)****: в JOURNAL.md, no notification to human unless signaled above
- ****End-of-experiment report****: day 21, full analysis, recommended next steps

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Scenario 1 — Positive progress:

[CDATA-A1] 14:30 | INFO
 $\alpha > 0$ confirmed: proliferating wells show 2.3x higher polyGlu vs contact-inhibited at day 14 (p=0.003, n=6 wells each). Keeping protocol. Full data: dashboard/cycle_156

Scenario 2 — Warning:

[CDATA-A1] 03:15 | WARN
Cell density dropped 28% in wells A4, A5 between 02:00 and 03:00 imaging. Possible contamination. Recommend visual check at 8 AM. Wells A1-A3, A6 normal. Continuing schedule, flagged for review.

Scenario 3 — Critical:

[CDATA-A1] 05:42 | CRIT
CO₂ dropped to 3.1% (target 5.0%). Opened valve 45s × 3, no recovery. Tank may be empty. HALT imaging. Cells survive ~4 hrs at low CO₂. Action needed before 10 AM. Checkpoint saved.

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Aim A.3 (CCP1/TLL6-OE causal): добавить signals for population doubling rates, transduction efficiency (GFP% cells), senescence onset timing.

Phase B (HSC mouse): полностью другой PROMPT.md — focus на chimerism, flow cytometry results, не microscopy.

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“Our automated microscopy platform pairs low-cost retrofit (\$4,500) with **prompt-driven AI supervision** (Claude Code /overnight). The experimenter articulates goals и thresholds в natural-language PROMPT.md; Claude continuously monitors data and signals the human only when strategically important events occur (positive progress, warnings, critical failures). This eliminates the ‘experimenter wake-up cost’ traditional time-lapse requires and democratizes AI-assisted experimental science for resource-limited labs. All PROMPT templates, tool functions, and decision logs released under MIT license concurrent with Phase A preprint.”

Follow-up upgrades (post-Phase A, if Go decision)

If Phase A passes и funding extends (Phase B or follow-on grant):

- **Hamamatsu ORCA-Fusion BT** upgrade (+\$4k used) → higher sensitivity for low-expression fluorescence
- **Perfect Focus System** (Nikon PFS) retrofit (+\$3k used) → hardware-based autofocus, no software drift

- **Motorized filter turret** (+\$1k) → automated channel switching, faster imaging
- **Liquid handling robot** (OpenTrons OT-2 \$2.5k new) → automated media change, ideal for 3-4 week contact-inhibition protocols

Budget for upgrades: place into Phase B or post-grant \$10-15k capital line.

Decision for LOI v24 submission (2026-04-24)

Recommendation: Вариант А (\$4,500 total equipment line) — максимально cost-effective, полностью покрывает Phase A experimental requirements (automated day-night imaging, environmental control, remote monitoring, data backup). Upgradeable to Phase B levels later.

Action:

1. Amend Phase A budget: “Microscope Upgrade” \$2k → “Microscope Automation & Upgrade” \$4.5k (reallocate \$2.5k from Consumables)
 2. Add this document as supplementary material
 3. Include in Cover Letter: “A custom-engineered automated time-lapse imaging platform (see Supplementary: Automated Microscopy Setup) will enable round-the-clock tracking of centriolar polyGlu dynamics in live fibroblasts, maximizing temporal resolution within budget constraints.”
-

Engineered 2026-04-21 for Impetus LOI v24 submission 2026-04-25.